

kinexus[®] Wi-Fi Gateway

Quick Start Guide

kinexus Wi-Fi Gateway Quick Start Guide

The kinexus Wi-Fi Gateway is used to transfer treatment data to your dialysis center. The data transfer only occurs when the Liberty Select Cyclor is connected to the Gateway and a Wi-Fi connection is established.

This guide is intended to provide patients with quick guidance to set up and connect the kinexus Wi-Fi Gateway to the Liberty Select Cyclor.

Please refer to the User Guide for warnings and more detailed instructions. This supplemental material does not replace the kinexus Wi-Fi Gateway User's Guide. Users should review and refer to the User's Guide for complete information on this product.

Have questions or issues? Contact Technical Support 800-227-2572



Before you begin:

- Turn on your mobile device and make sure Bluetooth is enabled
- Have the name of the intended Wi-Fi network available
- Have the Wi-Fi network password available

Setting up the Gateway

1. **Open** package and **check** for missing or damaged parts.



2. **Connect** the power cord to the Gateway and **plug** the power adapter into a wall outlet.




3. Wait approximately **5 minutes** for the **Status** LED to turn **green**. You may hear an audible signal “Gateway is ON”




Installing the kinexus application

1. Download and install the **kinexus Wi-Fi** application from the App Store or Google Play on your mobile device.


If your healthcare provider supplied a mobile device for your PD therapy, the app may already be installed. If so, skip to the next step.

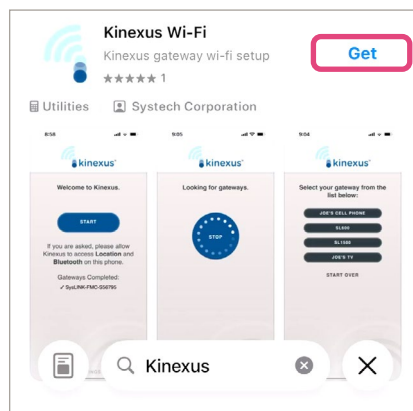


2. **Open** the app and allow kinexus access, if prompted.

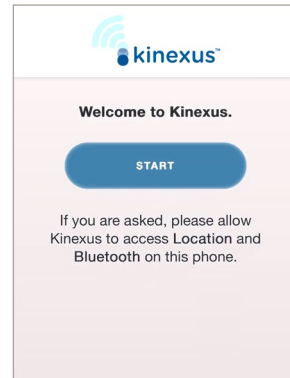


3. Tap **START**, ensure the Gateway is ready and then tap **CONTINUE**.

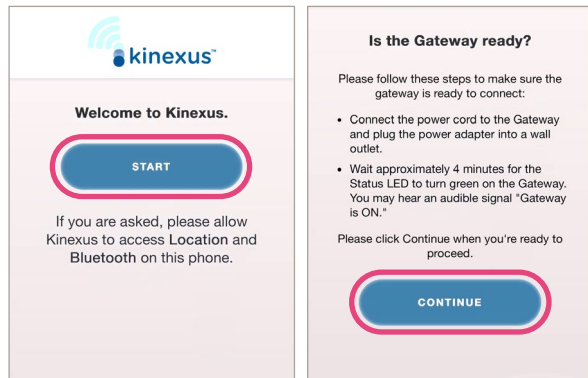


4. **Select** your Gateway and tap **Pair** when prompted

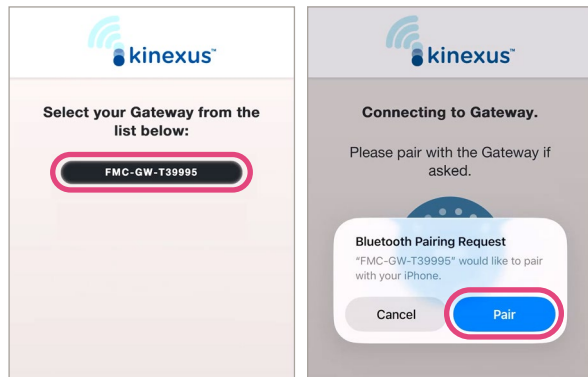

5. Select your home Wi-Fi network, **enter** your password then tap **APPLY** when prompted.

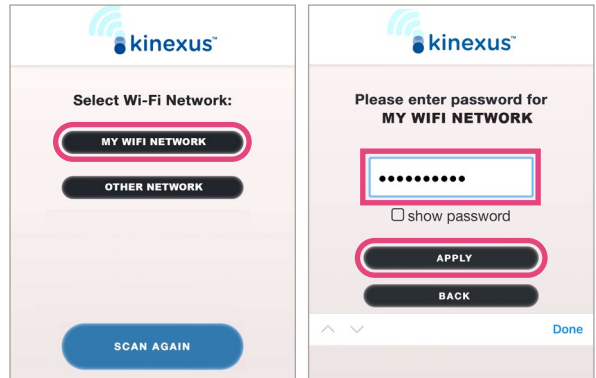

6. Success! appears when connection is complete.

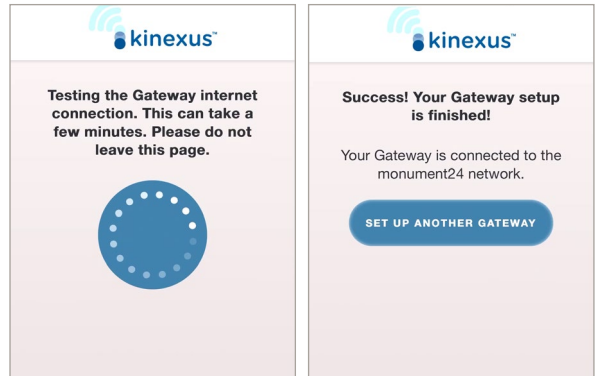

Connecting Liberty Select Cyclor to the Gateway

1. Ensure the Cyclor is powered **OFF**.



2. **Connect** the RS-232 end of the cable to the Cyclor. **Press** between the thumb screws to firmly seat the connector.



3. **Turn** both thumb screws clockwise until finger tight.




4. **Plug** the USB end of the cable into the USB port on the Gateway.

Note: The Wi-Fi Gateway should be powered ON for **5 minutes** before turning on the Cyclor.



5. **Turn on** the Liberty Select Cyclor.

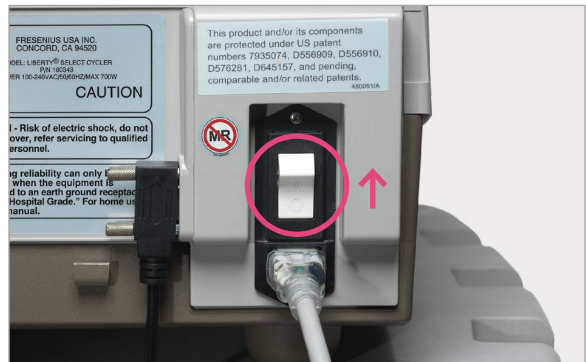



Read the Instructions for Use for safe and proper use of this device. For a complete description of hazards, contraindications, side effects, and precautions, see full package labeling at www.fmcna.com

Indications for Use: The Fresenius Liberty Select Cyclor is indicated for acute and chronic peritoneal dialysis.

Caution: Federal (US) law restricts this device to sale by or on the order of a physician.

Use: The kinexus Gateway is intended for transmission of home dialysis therapy data and patient data for patients using certain Fresenius Medical Care products, including peritoneal dialysis cyclor(s).

kinexus is not an electronic health record (EHR) or billing solution.

© 2026 Fresenius Medical Care. All Rights Reserved. Fresenius Medical Care, the triangle logo, Liberty, and kinexus are trademarks for which Fresenius Medical Care AG holds rights. All other trademarks are the property of their respective owners. P/N 480174 Rev C 01/2026